Ebola

Rapid Lateral Flow Test

Test rapide

Schnelltest

Rapid lateral flow test for the qualitative detection of the VP40 antigen of the Ebola virus in serum and throat swabs. For Professional In Vitro Diagnostic Use Only

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Ⓕ **Mode d'emploi**

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STADA Diagnostik Ebola Rapid Lateral Flow Test

Rapid lateral flow test for the qualitative detection of the VP40 antigen of the Ebola virus in serum and throat swabs

For Professional In Vitro Diagnostic Use Only

REF 2.13.1.122

Intended use

The STADA Diagnostik Ebola Rapid Lateral Flow Test is a rapid visually read chromatographic immunoassay for the qualitative detection of Ebola virus antigen in human serum and throat swabs.

It is used to detect Ebola infection in febrile conditions. It is intended to be used only together with clinical information/methods and other diagnostic in vitro procedures.

Summary

Ebola hemorrhagic fever (Ebola HF) is a severe, often fatal disease of humans and nonhuman primates (monkeys, gorillas and chimpanzees) that has appeared sporadically. The virus is one of two members of a family of RNA viruses called the Filoviridae. There are five identified subtypes of Ebola virus. 90 % of the outbreaks are caused by the Ebola Zaire strain.

Principle

The STADA Diagnostik Ebola Rapid Lateral Flow Test is a rapid, qualitative, membrane-based immunoassay for detection of the VP40 antigen of the Ebola virus in human serum and throat swabs.

After applying the sample to the sample well, it can react with particles coated with anti-Ebola antibodies. The mixture migrates upward on the membrane by capillary action.

Upon reaching the test line, antigen bound to the particles can react with the capture reagent, an anti-Ebola antibody immobilized on the test line of the membrane, and generate a colored line.

Appearance of a colored line in the test region indicates a positive result (presence of Ebola antigen), while its absence indicates a negative result.

For process control, a colored line will always appear in the control line region, indicating correct reaction conditions and procedure (proper volume of sample has been added and membrane wicking has occurred).

Contents of the kit

Kit for 20 determinations

20 test cassettes in aluminium pouch including desiccant bag

The kit comprises particles coated with anti-Ebola antibody, anti-Ebola antibodies and anti-mouse IgG immobilized on nitrocellulose membrane.

20 transfer pipettes

20 sterile polyester-tipped applicators/
swabs **STERILE EO** **CE** 0086,  Puritan Medical
Products Company LLC

20 extraction tubes

1 workstation for up to 8 extraction tubes

1 x 10 ml extraction buffer, ready-to-use 

The buffer solution contains detergent and preservative.

1 package insert

Equipment required but not provided

Timer, for measuring the time between application of the sample and interpretation of the results after ten minutes.

Sample collection device for serum.

Storage and stability

Store the STADA Diagnostik Ebola Rapid Lateral Flow Test unopened at +4 °C to +30 °C (+39 °F to +86 °F) and use it by the expiration date given on the label.

Do not use any component past the labeled expiration date.

Do not freeze.

The test cassette is sensitive to light and to humidity, therefore do not use test devices that have become wet or cassettes whose pouch is damaged or shows signs of damage!

After opening of the test cassette pouch, the test has to be performed immediately, within max. 1 hour.

After opening, the buffer vial can be stored at +4 °C to +30 °C (+39 °F to +86 °F) and safely be used for three months.

Samples

Samples should be collected, handled and shipped using appropriate infection control precautions for Ebola, other hemorrhagic fever viruses or other infectious diseases.

Shipping may be performed according to the policies of the sender, customs regulations, and the requirements of the receiving laboratory. If samples are to be shipped, they should be packed in compliance with statutory regulations for transportation of etiologic agents.

Serum

Only use human serum collected and prepared according to the instructions of the manufacturer of the sampling device for sample collection.

Only clear, non-hemolyzed samples may be used.

Serum may be stored at +2 °C to +8 °C (+35 °F to +46 °F) for up to 3 days. For long-term storage, samples may be kept below -20 °C/-4 °F.

Bring samples to room temperature prior to testing. Frozen samples must be completely thawed and mixed well prior to testing.

Samples should not be frozen and thawed repeatedly.

Throat swabs

Use sample at room temperature only (+15 °C to +30 °C/+59 °F to +86 °F).

Only use fresh throat swabs. Do not store. Do not freeze.

Warnings, notices and precautions

Handle all samples and used test components as though they contain infectious agents. Observe established precautions against microbiological hazards throughout the procedure. Appropriate biosafety practices must be complied with when handling samples and test components. These precautions include, but are not limited to, the following.

- Wear personal protective equipment: powder-free gloves, lab coat, eye protectors etc.
- Do not eat, drink, smoke, apply cosmetics, or handle contact lenses in areas where these materials are handled.
- Clean and disinfect all spills of samples using suitable disinfectant, such as 0.5 % sodium hypochlorite.
- Decontaminate and dispose of all samples, test components and other potentially contaminated materials in accordance with the local regulations.

Testing procedure

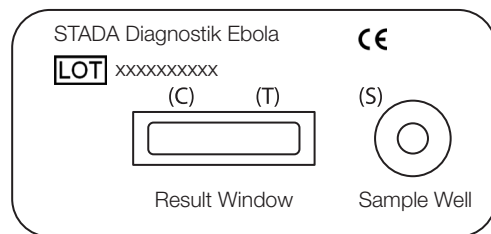
Fig. 1: Test cassette

Preparation

Allow the kit components and sample to equilibrate to room temperature (+15 °C to +30 °C/+59 °F to +86 °F) prior to testing.

- Remove the test cassette from the sealed pouch and use it as soon as possible. The test should be performed within 1 hour after opening the pouch.

- Place the test cassette on a clean and level surface.
- Clearly indicate the patient/sample identification on the test device to avoid confusion.
- Place the extraction tube in the workstation and indicate the patient/sample identification on the tube.

Throat swabs

- Hold the extraction buffer bottle vertically and add **8 to 10 drops of the extraction buffer solution to the extraction tube.** Only use the buffer solution provided with this test kit.
- Insert the provided sterile swab into the patient's oral cavity.
- Very gently push the swab until resistance is met at level of the pharynx.
- Gently rotate the swab for a few times over the mucous membrane.
- Immediately insert the swab with collected throat swab sample into the tube, compress the bottom of the tube and rotate the swab vigorously to mix the reagents for at least 1 minute. Press the swab against the inside of the tube, and squeeze as much liquid as possible from the swab by pressing and rotating the fiber portion against the wall of the tube. Discard the swab.
- Let the mixture incubate for 1 to 15 minutes.
- Mix the contents of the tube again by gentle swirling.
- Now the mixture is ready to be tested. Fill the provided transfer pipette with the mixture.
- Hold the provided transfer pipette over the sample well of the test cassette and **squeeze 6 drops into the well.** Wait until each drop is absorbed before adding additional drops.
- Start the timer.
- Wait 10 minutes from addition of the sample/buffer mixture and then read the results. **Warning:** Do not interpret the results after more than 15 minutes, to avoid misinterpretation.

Serum

- Hold the extraction buffer bottle vertically and add **3 drops (approx. 100 µl) of the extraction buffer to the extraction tube.** Only use the buffer included in the test kit.
- Hold the provided transfer pipette vertically, transfer **3 drops of sample (approx. 100 µl) into the extraction tube** and incubate for 1 to 15 minutes.
- Now the mixture is ready to be tested. Fill the provided transfer pipette with the mixture, hold the pipette vertically over the sample well of the test cassette and drop the whole mixture into the sample well. Wait until each drop is absorbed before adding additional drops.
- Start the timer.
- Wait 10 minutes from addition of the sample/buffer mixture and then read the results. **Warning:** Do not interpret the results after more than 15 minutes, to avoid misinterpretation.

Quality control

A process control is included in the test. A colored line appearing in the control region (C) is the internal process control. It confirms sufficient sample volume and correct technique.

If the control line fails to appear (insufficient sample volume or incorrect techniques are the most likely reasons), the test result is INVALID (figure 4). Check the procedure and repeat the test with a new test cassette. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

Interpretation of results

POSITIVE: Two distinct colored lines appear. One line should be in the control region (C), and another line should be in the test region (T) (figure 2).

Any visible colored line in the test region (T) should be interpreted as positive, even if the line in the test region (T) appears lighter or darker than the line in the control region (C).

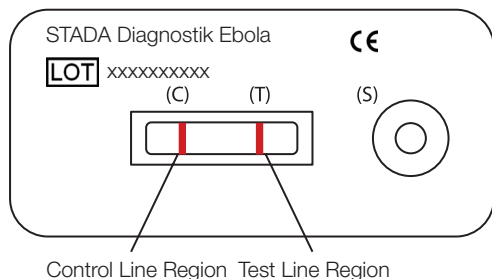
Fig. 2: Schematic illustration of a positive result

NEGATIVE: One colored line appears in the control region (C). No visible line appears in the test region (T) (figure 3).

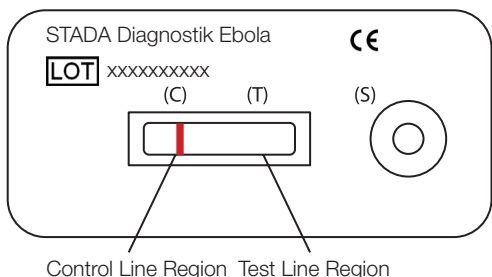
Fig. 3: Schematic illustration of a negative result

INVALID: Control line fails to appear (figure 4). If there is no colored line in the control region (C), even if a colored line appears in the test region (T), the result is invalid, and the test should be repeated with a new test cassette.

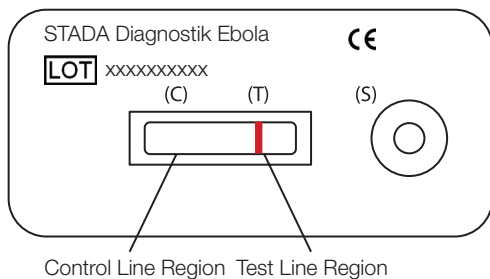
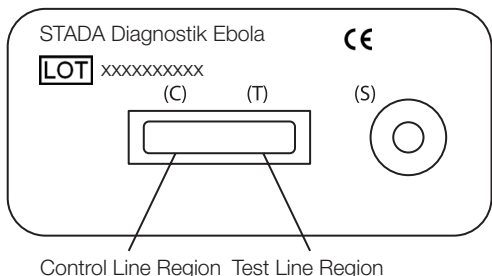
Fig. 4: Schematic illustration of invalid results

Limitations

The STADA Diagnostik Ebola Rapid Lateral Flow Test is designed for qualitative in vitro detection of the VP40 antigen of the Ebola virus in human serum and throat swabs. Other body fluids and pooled samples may not give accurate results, and should not be used.

The STADA Diagnostik Ebola Rapid Lateral Flow Test should be used only on samples from patients with suspected Ebola infections (febrile diseases) as an aid to detect Ebola hemorrhagic fever. It should not be used as the sole criterion in diagnosing Ebola infections.

The STADA Diagnostik Ebola Rapid Lateral Flow Test will indicate presence of Ebola virions at densities of $\geq 10^6$ PFU/ml in the sample. Therefore it should not be used to exclude Ebola infections.

Results should be interpreted by a physician only.

Results should be used only in synopsis with clinical information/symptoms/methods and other In Vitro Diagnostic procedures (e.g. PCR techniques).

The intensity of the colored line in the test region does not correlate to the titer of antigen in the sample.

A negative result does not exclude the possibility of Ebola infection at any time.

A positive result should be confirmed using a different method, and the result should be evaluated together with the overall clinical evaluation.

False-positive or false-negative results may be produced by:

- insufficient sample volume or incorrect techniques
- use of inappropriate sample or materials (other buffer than the one included in the test kit, other volumes of samples and buffer than those indicated)
- disregard of the incubation time of 1 minute of sample and extraction buffer
- disregard of the waiting time of 10 minutes after applying the mixture to the test cassette

Low titers of Ebola may result in a faint line appearing in the test region after a prolonged time.

The "Hook effect" in case of very high levels of antigen in the sample can lead to false-negative results.

Apart from EDTA, no substance has been demonstrated to trigger cross-reactions.

Performance characteristics

The performance of STADA Diagnostik Ebola Rapid Lateral Flow Test has been compared with that of the PCR system RealStar® Ebola virus RT-PCR Kit 1.0 based on samples of African patients from Conakry in Guinea.

		STADA Diagnostik Ebola		total
		+	-	
PCR	+	12	1	13
	-	2	103	105
total		14	104	

Analytic sensitivity: $\leq$ CT 23

Specificity: 98 %

References

1. Lucht A, Grunow R, Möller P, Feldmann H, Becker S: Development, characterization and use of monoclonal antibodies for the detection of Ebola virus. *J Virol Methods* 2003; 111:21-8
2. Grolla A, Lucht A, Dick D, Strong JE, Feldmann H: Laboratory diagnosis of Ebola and Marburg hemorrhagic fever. *Bull Soc Pathol Exot* 2005; 98:205-9
3. Becker S, Feldmann H, Will C, Slenczka W: Evidence of occurrence of filovirus antibodies in humans and imported monkeys: Do sub-clinical filovirus infections occur worldwide? *Med Microbiol Immunol* 1992; 181:43-55
4. World Health Organization. Outbreak(s) of Ebola haemorrhagic fever in the Republic of Congo, January-April 2003. *Wkly Epidemiol Rec* 2003; 78:285-9
5. Hartmann H, Lübbers B, Cararetto M, Bausch W, Klos A, Köhl J. Rapid quantification of C3a and C5a using a combination of chromatographic and immunoassay procedures. *J Immunol Methods* 1993; 166:35-44
6. Ksiazek TG, Rollin PE, Jahrling PB, Johnson E, Dalgard DW, Peters CJ: Enzyme immuno-sorbent assay for Ebola virus antigens in tissues of infected primates. *J Clin Microbiol* 1992; 30:947-50
7. National Committee for Clinical Standards Clinical Waste Management: Approved Guideline. NCCLS Document GPS-A. Villanova, PA: NCCLS 1993; 13:1-18, 29-42
8. Towner JS, Rollin PE, Bausch DG, Sanchez A, Crary SM, Vincent M, Lee WF, Spiropoulou CF, Ksiazek TG, Lukwiya M, Kaducu F, Downing R, Nichol ST: Rapid diagnosis of ebola hemorrhagic fever by reverse transcription-PCR in an outbreak setting and assessment of patient viral load as a predictor of outcome. *J Virol* 2004; 78:4330-41

STERILE EO

Package contains a sterile product

CE 0086,  Puritan Medical Products Company LLC

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Références

1. Lucht A, Grunow R, Möller P, Feldmann H, Becker S: Development, characterization and use of monoclonal antibodies for the detection of Ebola virus. J Virol Methods 2003; 111:21-8
2. Grolla A, Lucht A, Dick D, Strong JE, Feldmann H: Laboratory diagnosis of Ebola and Marburg hemorrhagic fever. Bull Soc Pathol Exot 2005; 98:205-9
3. Becker S, Feldmann H, Will C, Slenczka W: Evidence of occurrence of filovirus antibodies

- in humans and imported monkeys: Do sub-clinical filovirus infections occur worldwide? Med Microbiol Immunol 1992; 181:43-55
4. World Health Organization. Outbreak(s) of Ebola haemorrhagic fever in the Republic of Congo, January-April 2003. Wkly Epidemiol Rec 2003; 78:285-9
5. Hartmann H Lübbers B, Cararetto M, Bautsch W, Klos A, Köhl J. Rapid quantification of C3a and C5a using a combination of chromatographic and immunoassay procedures. J Immunol Methods 1993; 166:35-44
6. Ksiazek TG, Rollin PE, Jahrling PB, Johnson E, Dalgard DW, Peters CJ: Enzyme immuno-sorbent assay for Ebola virus antigens in tissues of infected primates. J Clin Microbiol 1992; 30:947-50
7. National Committee for Clinical Standards Clinical Waste Management: Approved Guideline. NCCLS Document GPS-A. Villanova, PA: NCCLS 1993; 13:1-18, 29-42
8. Towner JS, Rollin PE, Bausch DG, Sanchez A, Crary SM, Vincent M, Lee WF, Spiropoulou CF, Ksiazek TG, Lukwiya M, Kaducu F, Downing R, Nichol ST: Rapid diagnosis of ebola hemorrhagic fever by reverse transcription-PCR in an outbreak setting and assessment of patient viral load as a predictor of outcome. J Virol 2004; 78:4330-41

STERILE EO

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CE 0086, Puritan Medical Products Company LLC

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Analytische Sensitivität: $\leq CT 23$

Spezifität: 98 %

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1. Lucht A, Grunow R, Möller P, Feldmann H, Becker S: Development, characterization and use of monoclonal antibodies for the detection of Ebola virus. J Virol Methods 2003; 111:21-8
2. Grolla A, Lucht A, Dick D, Strong JE, Feldmann H: Laboratory diagnosis of Ebola and Marburg hemorrhagic fever. Bull Soc Pathol Exot 2005; 98:205-9
3. Becker S, Feldmann H, Will C, Slenczka W: Evidence of occurrence of filovirus antibodies in humans and imported monkeys: Do sub-clinical filovirus infections occur worldwide? Med Microbiol Immunol 1992; 181:43-55
4. World Health Organization. Outbreak(s) of Ebola haemorrhagic fever in the Republic of Congo, January-April 2003. Wkly Epidemiol Rec 2003; 78:285-9
5. Hartmann H, Lübbers B, Cararetto M, Bautsch W, Klos A, Köhl J. Rapid quantification of C3a and C5a using a combination of chromatographic and immunoassay procedures. J Immunol Methods 1993; 166:35-44
6. Ksiazek TG, Rollin PE, Jahrling PB, Johnson E, Dalgard DW, Peters CJ: Enzyme immuno-sorbent assay for Ebola virus antigens in tissues of infected primates. J Clin Microbiol 1992; 30:947-50

7. National Committee for Clinical Standards
Clinical Waste Management: Approved
Guideline. NCCLS Document GPS-A. Villano-
va, PA: NCCLS 1993; 13:1-18, 29-42
8. Towner JS, Rollin PE, Bausch DG, Sanchez
A, Crary SM, Vincent M, Lee WF, Spiropou-
lou CF, Ksiazek TG, Lukwiya M, Kaducu F,
Downing R, Nichol ST: Rapid diagnosis of
ebola hemorrhagic fever by reverse transcrip-
tion-PCR in an outbreak setting and assess-
ment of patient viral load as a predictor of out-
come. J Virol 2004; 78:4330-41

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Distributed by/Distribué par/Vertrieb durch: STADApHarm GmbH, 61118 Bad Vilbel, Germany, www.stada.de